



User Manual





Table of contents

01 What is billy?

02 Product overview

03 How to set it up

04 APP features

05 Folding instructions

06 Warranty

07 FAQ





What is illy?

"Healthy Kids, Happy Parents – Your Smart Medication Buddy"


The Billy is a smart Medication Dispenser that helps children and parents making medicine intake easier.

By using the Billy Parents don't have to worry about supervising their kids while taking their medication, the included camera system makes controll easier and our integrated App system even creates the option for supervision from work or somewhere else.


Moreover children won't keep fighting against taking medication, the design and integrated rewarding system will make it as appealing as playing a game!



Product Overview

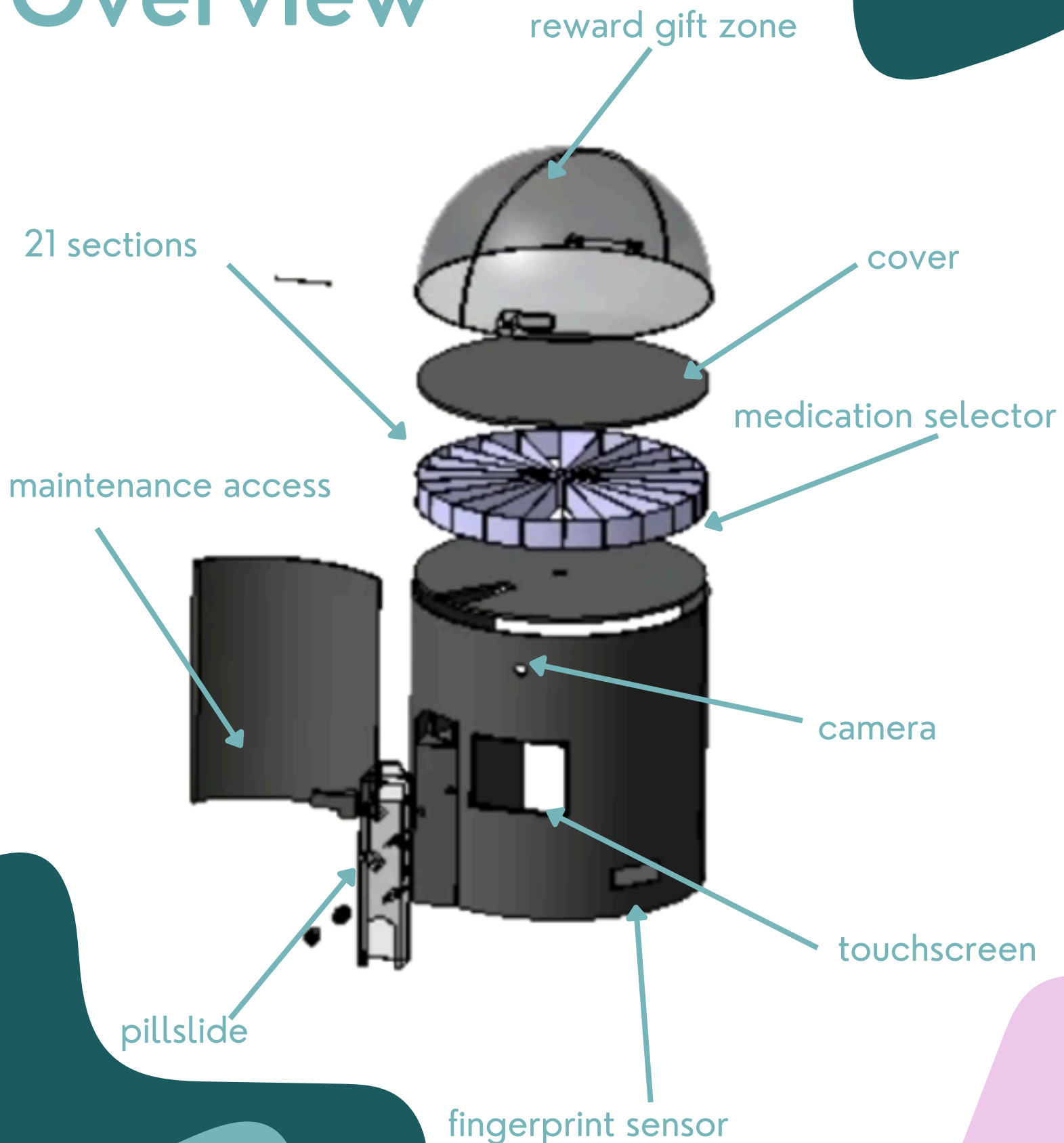

How to set it up

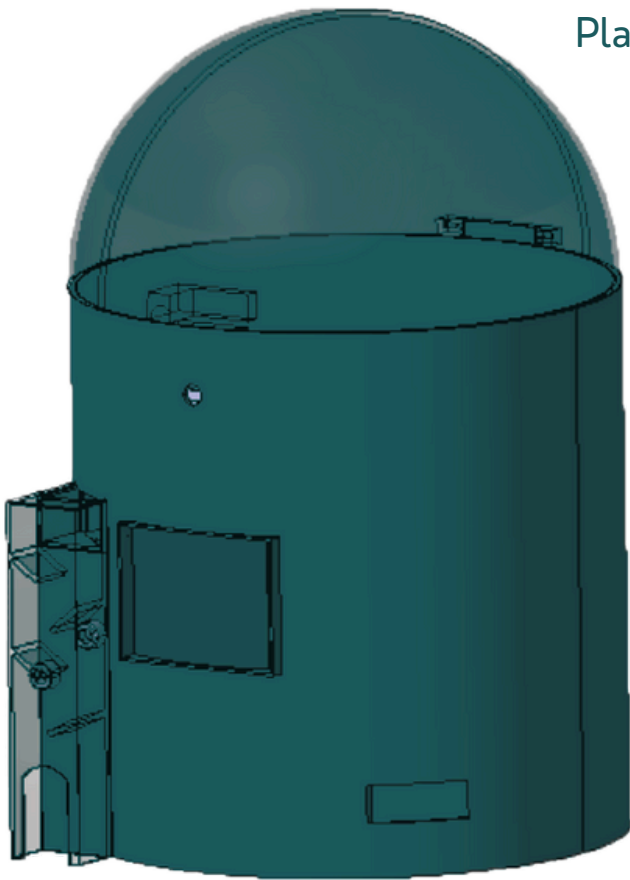

Place it on a steady surface

Plug it in

For using the Billy WI-FI access is necessary

Connect to the app with the following instructions


Fill up the rewarding area

Possible rewards:

- put a activity on a paper and into a ball to make it seem like a candy machine, vouchers
- sweets are always an option

Connect to the App

How to insert the medication

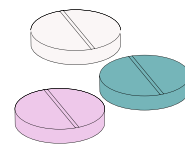

Parent fingerprint authentication to get access


Parentsmode is activated

Access to the pill mechanism is possible

Put in the medication

Close the door

Put your fingerprint again to close it up

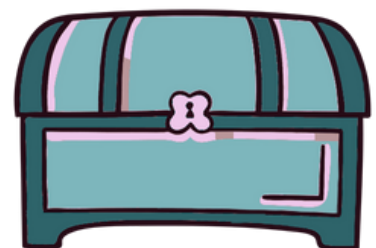


Give the packaging a second life
Instructions on page 8

APP Features




Web APP available

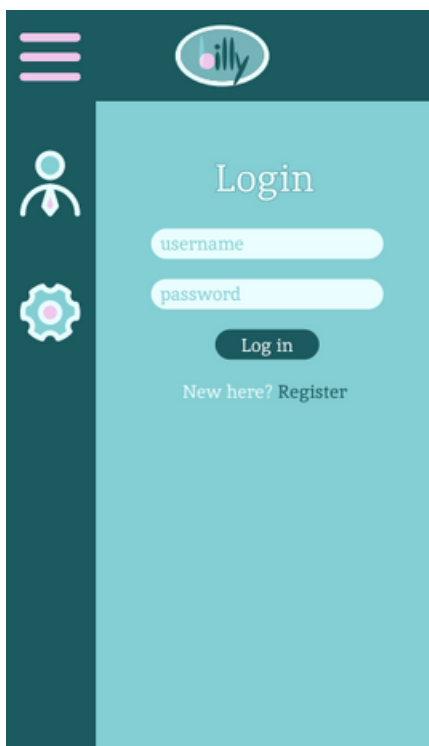

1. Register and create your account



2. Main Menu

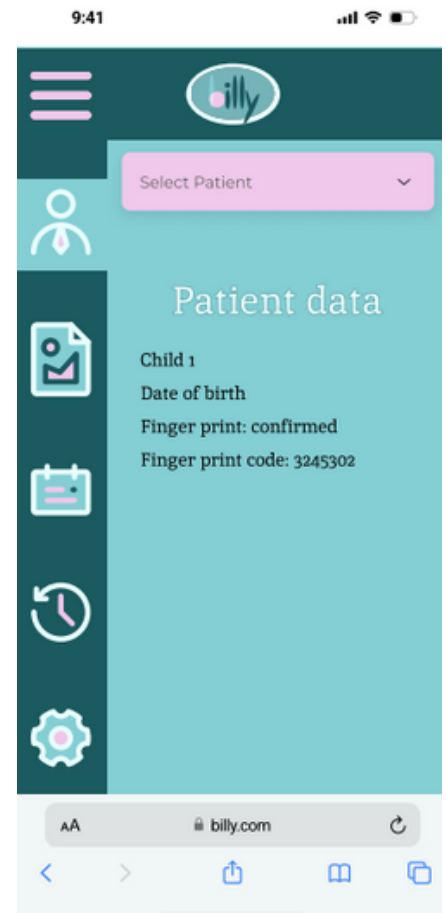


3. Patient data, to see the collected Data of chosen client


APP Features



4. Photo Library, see taken pictures

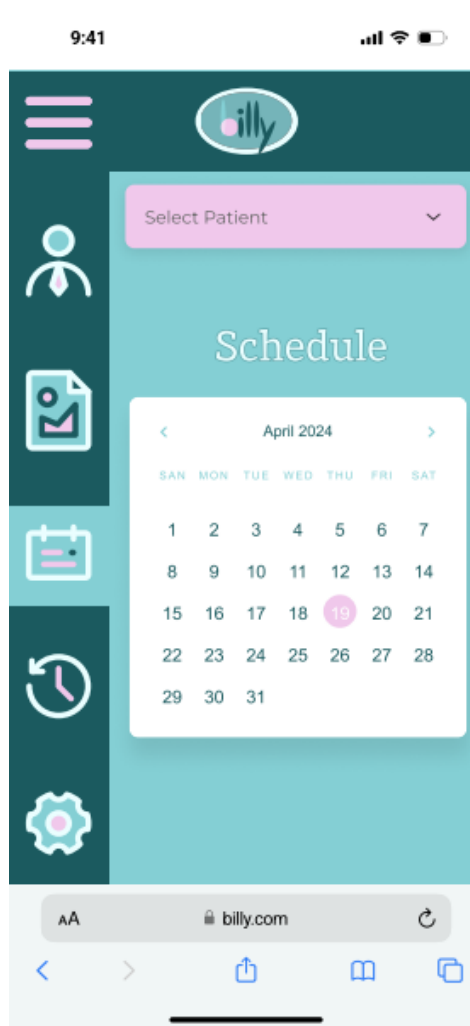

5. Calendar, choose patient, set and control the created shedule. Also changes can be made here

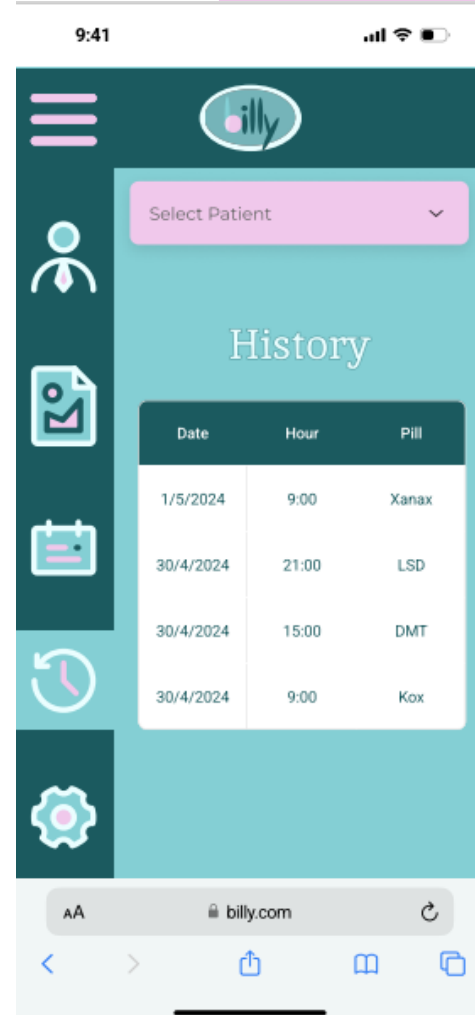



6. History, choose patient and see the history of medication intake

Scan the code to get the direct link to the web application



Folding Instructions

Give the Packaging a new life
by creating your own treasure chest

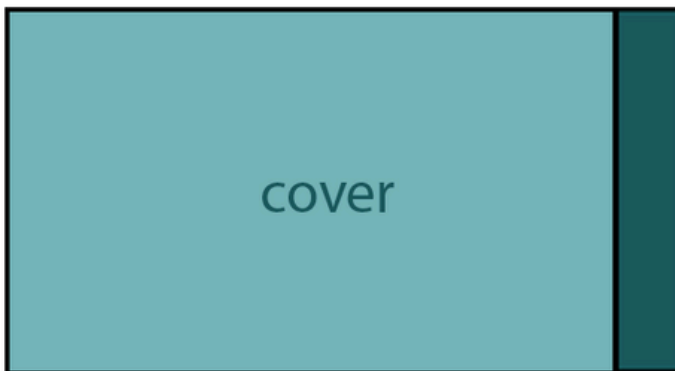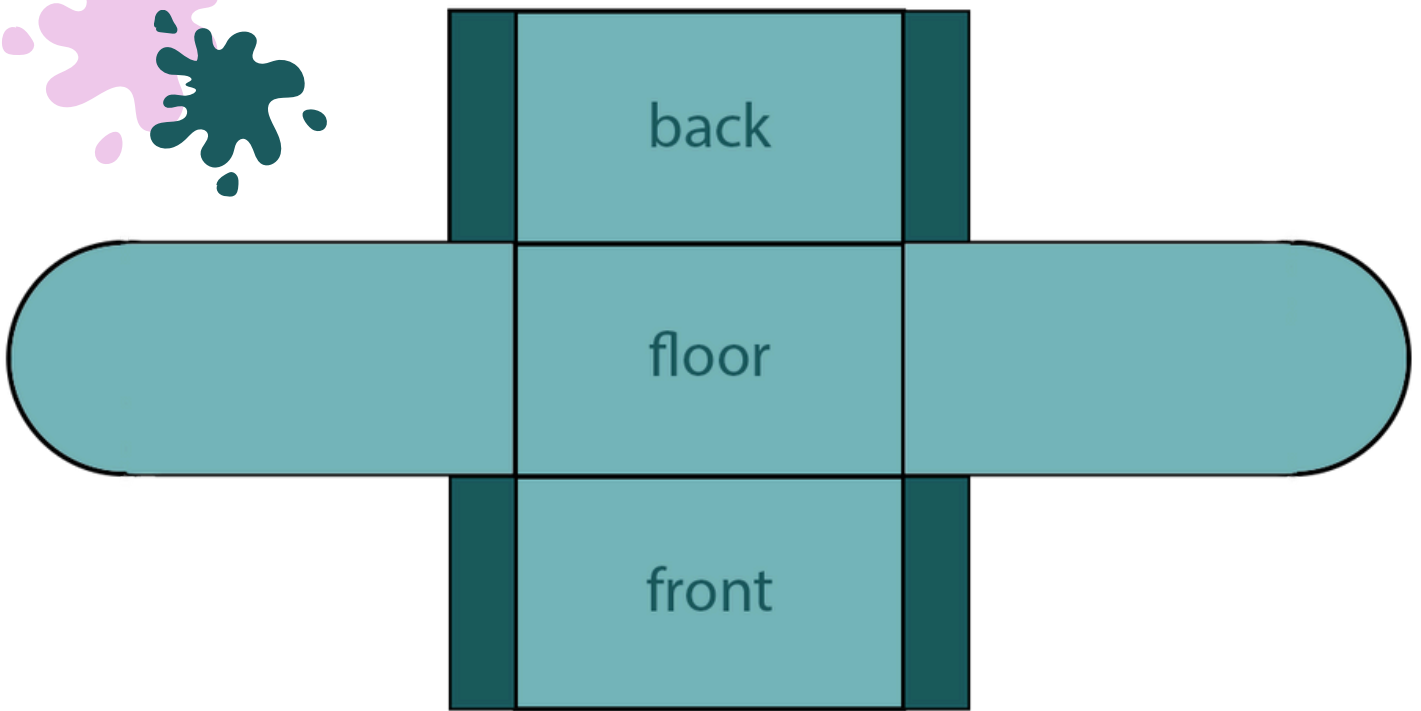

1. cut out the given guiding lines
2. align the sides
3. glue the ends together



and now it's time for the fun part
get creative and design your own
masterpiece



Frequently Asked Questions

01

Renting

Require the Billy for a restricted duration? Rest assured, we offer rental options for the Billy as well. Visit our website for further details.

02

Disposal

As the Billy contains various materials, we provide the option to return it or deliver it to one of our authorized retail outlets, ensuring proper disposal.

03

Technology

In the cause of technical issues with the device or challenges with the user interface of the web application, we kindly encourage you to contact our customer service hotline or refer to our website for assistance.


www.billy.com/interface-help

04

Broken Parts

Should any component of the Billy become damaged, please do not hesitate to contact our customer service team. We maintain a supply of replacement parts to address such occurrences promptly.



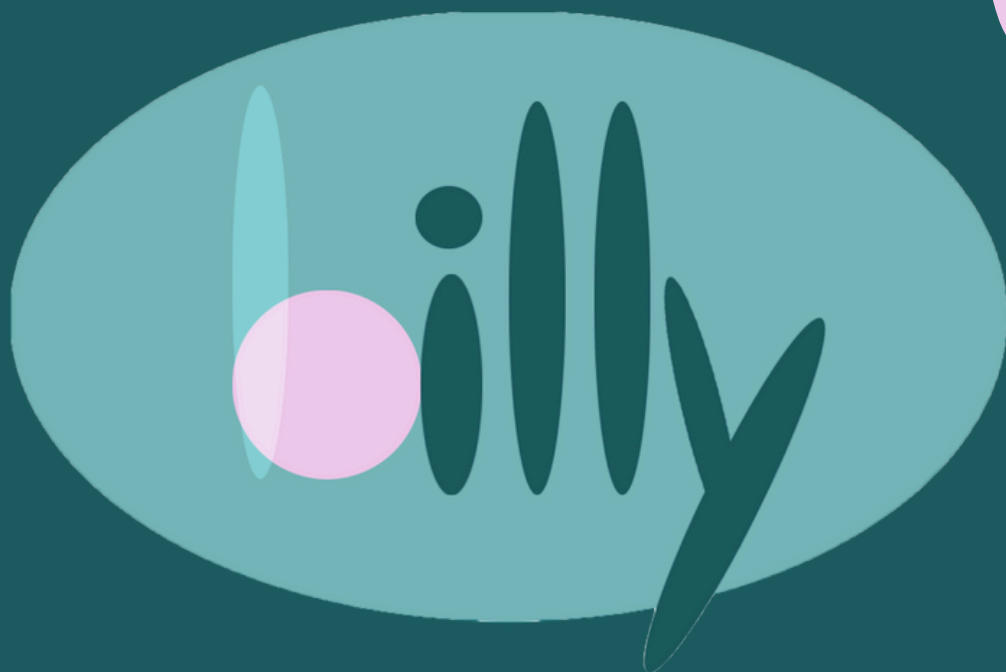

Contact us!



Porto, Portugal



+351 12322133



www.EPS2024-billy.com



Łukasz Borowski, Tibo Clauwaert, Lena Ehrenhofer
Tamara Kronshagen, Noé Oliveira, Stijn Steyaert

